




Organ Recovery systems

LifePort® Kidney Transporter Disposable Oxygenation Perfusion Circuit

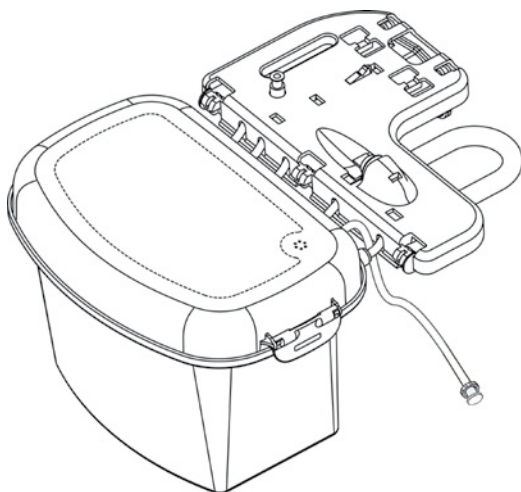

Instructions for Use

LifePort Kidney Transporter Disposable Oxygenation Perfusion Circuit is for use only with LifePort Kidney Transporter system.

 **Organ Recovery Systems, Inc.**
One Pierce Place, Suite 475W
Itasca, IL 60143
USA

T +1.847.824.2600
F +1.847.824.0234
Perfusion Helpline
+1.866.682.4800

Organ Recovery Systems NV
Culliganlaan 1B
1831 Diegem
Belgium

T +32.2.715.0000
F +32.2.715.0009
Perfusion Helpline
+32.2.715.0005

ORS Representacoes do Brasil Ltda.
170 Moema Avenue, Suite 11 & 12
Sao Paulo, SP 04077-020
Brazil

T +55.11.9.8638.0086

Perfusion Helpline
+55.11.9.8638.0086

LKT201X

INDICATIONS FOR USE

LifePort Kidney Transporter (LKT) system is intended to be used for the continuous hypothermic machine perfusion of kidneys for the preservation, transportation and eventual transplantation into a recipient.

DEVICE DESCRIPTION


LifePort Kidney Transporter Disposable Oxygenation Perfusion Circuit is used to contain the kidney and perfusate under aseptic conditions during transport, with the option to introduce oxygen into the perfusate.

The complete LifePort Kidney Transporter system is comprised of the following components:


- LifePort Kidney Transporter (LKT100P/LKT101P/LKT101PNG)
- LifePort Kidney Transporter Disposable Oxygenation Perfusion Circuit (LKT201X)
 - Oxygenation Delivery Tube (900-00044)
 - Oxygenation Connector (900-00045)
- LifePort Kidney Transporter Disposable Sterile Drape (LKT300)
- LifePort Kidney Transporter Disposable Cannula (CAN/UCAN)

INSTRUCTIONS FOR USE


Prepare LifePort Kidney Transporter as described in LifePort Kidney Transporter Operator's Manual. Prepare KPS-1®, Kidney Perfusion Solution in accordance to the Instructions for Use. Maintain perfusate at proper temperature, approximately 2–8°C.

 **WARNING:** Where noted, perform the following procedure on an aseptic field using aseptic technique.


1. *Using standard aseptic technique*, prepare a sterile field on a work table and introduce all necessary materials.
2. *Using standard aseptic technique*, remove LifePort Kidney Transporter Disposable Perfusion Circuit, Oxygenation Delivery Tube, and Oxygenation Connector from their respective packaging.

 **WARNING:** Visually check LifePort Kidney Transporter Disposable Products before use. Do not use if parts are cracked, broken, or disconnected.


3. *Using standard aseptic technique*, remove Outer Perfusion Circuit Lid and Inner Perfusion Circuit Lid and place onto the sterile field.
4. *Using standard aseptic technique*, remove the Kidney Cradle and set aside within the sterile field.

 **CAUTION:** Verify the Uptake Line is within the Low Fluid Detection Chamber inside the Organ Cassette.


5. *Using standard aseptic technique*, attach the Oxygenation Delivery Tube and Oxygenation Connector to Oxygenation Perfusion Circuit.

 **WARNING:** The Oxygenation Connector must be applied to the Oxygen Delivery Tube to prevent leakage or contamination of the aseptic field.


6. *Using standard aseptic technique*, fill Oxygenation Perfusion Circuit with 1 Liter chilled (2–8°C) perfusate.
7. *Using standard aseptic technique*, replace and secure Inner Perfusion Circuit Lid followed by the Outer Perfusion Circuit Lid.

 **WARNING:** Oxygenation Perfusion Circuit inner surfaces are considered sterile, while outer surfaces are not considered sterile.

8. Place Oxygenation Perfusion Circuit into LifePort Kidney Transporter.
9. Position the Tubeframe upright, perpendicular to the Pump Deck. Insert the hinges into the receivers before rotating flat onto the Pump Deck.

 **CAUTION:** Make sure Oxygenation Delivery Tube is not pinched or below the Tubeframe and is easily accessed from the front.

10. Open the Pumphead Raceway and stretch the Pump Tubing Loop around the Infusion Pump. Close and latch the Pumphead Raceway.
11. Rotate the Pump Deck Locking Arm 90 degrees until it clicks into place.
12. Connect the Pressure Sensor Cable from the Pump Deck to the Pressure Sensor Connector on the Tubeframe.
13. Remove and set aside the Oxygenation Connector from the Oxygenation Delivery Tube.
14. Connect oxygen supply tubing (not included) from an oxygen source with a regulator to the Oxygenation Delivery Tube.

 **CAUTION:** Only use medical-grade oxygen that meets the regulatory requirements for compressed medical gas.

NOTE: Oxygen supply tubing is not provided by Organ Recovery Systems.

15. Press and hold the **POWER** button until you hear an audible beep, then release.
 16. Press the **WASH** button to begin **WASH** mode.
 17. Turn on the oxygen and administer at a rate of 0.5 Liters/min for a minimum of 20 minutes.
-



CAUTION: During perfusate oxygenation, do not exceed an oxygen flow of 0.5 Liters/min to avoid overpressurizing the Organ Chamber while both Lids are secured.



CONSULT INSTRUCTIONS FOR USE: Follow the procedure outlined in the LifePort Kidney Transporter Disposable Cannula Instructions for Use (IFU) to cannulate and secure the kidney.

18. Once the kidney is ready to be placed into the Organ Chamber, turn off oxygen and disconnect the oxygen supply tubing from the Oxygenation Delivery Tube and oxygen source.
 19. Disconnect the Oxygenation Delivery Tube from the Oxygenation Perfusion Circuit.
 20. Fasten the Oxygenation Connector to the open Oxygen Delivery Tube on the Oxygenation Perfusion Circuit to maintain aseptic conditions.
-



WARNING: The Oxygenation Connector must be applied to the Oxygen Delivery Tube to prevent leakage or contamination of the aseptic field.

21. Remove the Outer Perfusion Circuit Lid.
-



CONSULT INSTRUCTIONS FOR USE: Refer to LifePort Kidney Transporter Disposable Sterile Drape IFU to maintain aseptic conditions.

22. *Using standard aseptic technique*, remove Inner Perfusion Circuit Lid.
 23. *Using standard aseptic technique*, transfer the cannulated kidney in the Kidney Cradle to LifePort Kidney Transporter, being careful to avoid catching the Infuse Line.
-



CONSULT INSTRUCTIONS FOR USE: Refer to LifePort Kidney Transporter Operator's Manual to start the perfusion process and check for leaks.

24. *Using standard aseptic technique*, replace and secure Inner Perfusion Circuit Lid.
-



CONSULT INSTRUCTIONS FOR USE: Follow the procedure outlined in the LifePort Kidney Transporter Disposable Sterile Drape IFU for removal from LifePort Kidney Transporter.

25. Replace and secure Outer Perfusion Circuit Lid.
-



CONSULT INSTRUCTIONS FOR USE: Follow the procedure in LifePort Kidney Transporter Operator's Manual to continue the perfusion process.

STORAGE CONDITIONS

Store between 2°C and 40°C. Avoid excessive heat and humidity. Keep dry and keep away from direct sunlight. Sterile unless package is damaged or open.

CONTRAINDICATIONS

There are no known contraindications when used as directed.

TECHNICAL ASSISTANCE

Please contact Organ Recovery Systems 24/7 Perfusion Helpline at +1.866.682.4800 (toll free in US), +32.2.715.0005 (Belgium), +55.11.98638.0086 (Brazil), or +33.967.23.00.16 (France).

WARNINGS AND PRECAUTIONS



CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician.













CAUTION: LifePort Kidney Transporter Disposable Products should be stored indoors in a dry location out of direct sunlight.

















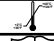

WARNING: Visually check LifePort Kidney Transporter Disposable Products. Do not use if parts are cracked, broken, or disconnected.



WARNING: For single use only. Do not reuse, reprocess or resterilize. Reusing, reprocessing, or resterilization of single-use devices creates a potential risk of patient or user infections due to contamination. This contamination may lead to injury, illness, or other serious patient complications.

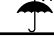-  **WARNING:** Use standard aseptic technique and universal precautions (e.g., gloves, masks, gowns, goggles or equivalent eye protection, biohazard bags) when handling the kidney, and when handling and disposing of LifePort Kidney Transporter Disposable Products and perfusate to prevent the possible transmission of pathogens to medical personnel and patients. The single practitioner, working alone must pay special attention to maintain these conditions.
-  **WARNING:** Before the operator begins to apply the LifePort Kidney Transporter Disposable Sterile Drape, the Outer Perfusion Circuit Lid should be removed. The exterior surfaces of the Outer Perfusion Circuit Lid are not considered part of a sterile field. Use standard aseptic technique when handling the interior surface.
-  **WARNING:** Oxygen must be handled in accordance with institutional safety procedures.
-  **WARNING:** Only experienced and properly instructed persons should handle oxygen.
-  **WARNING:** A fire may occur if oxygen comes in contact with a spark or other material that burns easily.
-  **WARNING:** Do not use electrical equipment that can produce a spark as it may ignite and burn when oxygen is in use.
-  **WARNING:** Do not smoke, light matches or use a lighter in a room where oxygen is being used.
-  **WARNING:** Do not use petroleum-based products.
-  **WARNING:** Use only oxygen approved lubricants and oxygen approved sealings.
-  **WARNING:** Room should be labeled with "Oxygen in Use" sign.

EXPLANATION OF SYMBOLS

	Warning/Caution		Use By, YYYY-MM-DD		Temperature Limits
	Lot Number		Date of Manufacture, YYYY-MM-DD		Consult Instruction for Use
	Reference Number		Manufacturer		Keep Away From Sunlight
	Do Not Reuse		Do Not Resterilize		Keep Dry
	Sterile Medical Devices Using Ethylene Oxide		Oxidizing Material		Prescription Medical Device
	Medical Device				



European Authorized Representative
Emergo Europe
Prinsessegracht 20
2514 AP The Hague
The Netherlands